


# HTTP-Request

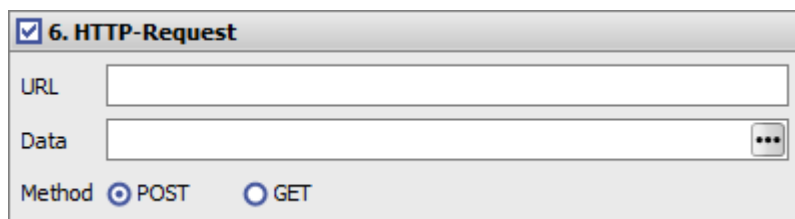

## Parameter URL

In field URL, an HTTP address can be specified. This will be called when the tool is started or the event occurs. The use of wildcards is possible.

```
http://LocalHost:9000/SETPAT -Client=%F -SelectImage="%7"
```

The following HTTP request is sent:

```
http://LocalHost:9000/SETPAT -Client=15 -SelectImage="10000R6I6"
```

The receiver must accept and interpret this call.

## Parameter Data

A text file can be configured in the Data field. This will be transferred with the HTTP request. The use of placeholders is also possible here. For example, a JASON or XML file can be filled with patient data and then sent.

Example: With the following JASON file template:

```
{
  "patientId" : "%0" ,
  "patientFirstName" : "%3"
}
```

wurde olgender Text übertragen:

```
{
  "patientId" : "10000" ,
  "patientFirstName" : "Petra"
}
```

The receiver must accept and interpret this call.

From:  
<https://www.onyxwiki.net/> - [OnyxCeph<sup>3</sup>™ Wiki]

Permanent link:  
<https://www.onyxwiki.net/doku.php?id=en:6httprequest>

Last update: **2023/10/09 09:01**

