

Module Add Image 3D



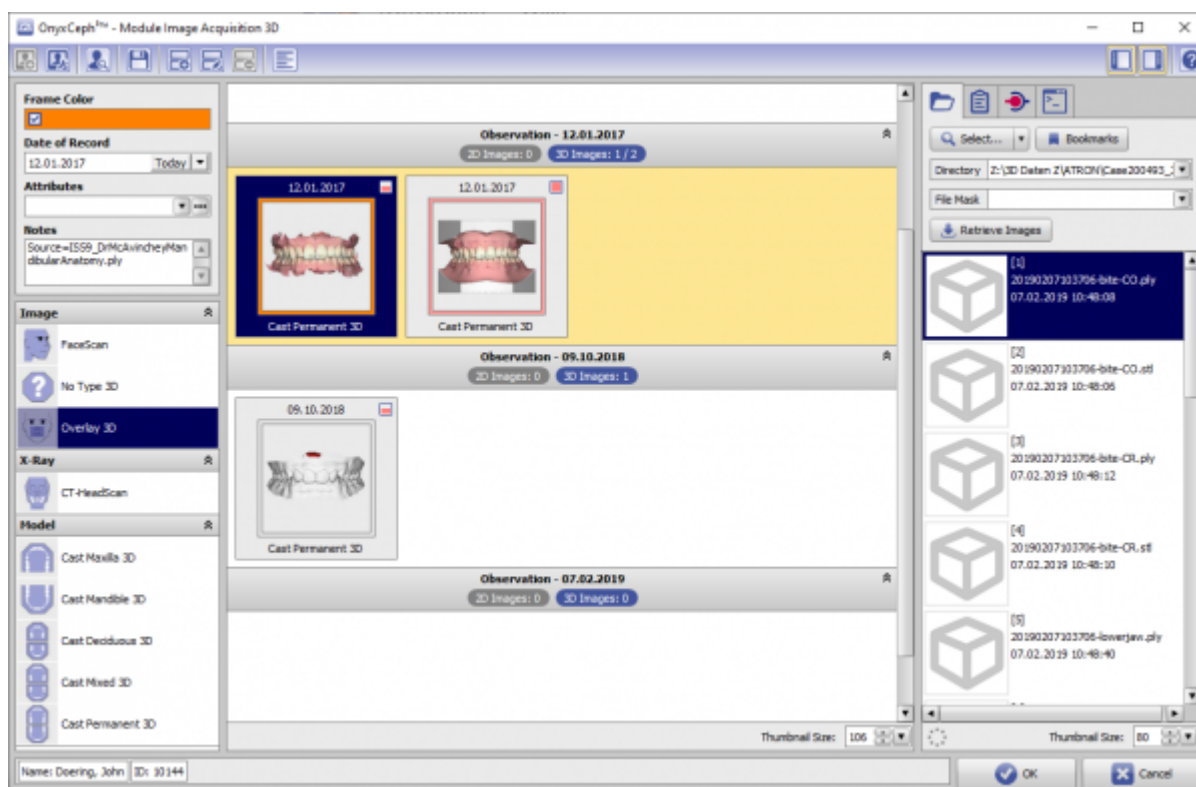
Program version	BASIC	2D PRO	3D PRO	LAB	OMS
Availability		✗			

In Module Add Image 3D different types of 3D data sets can be imported from pre-defined image sources, classified acc. to [session](#) and [image type](#) and completed with [attributes and comments](#) before saved to the Onyx database.





First Steps

1. Activate "Clipping (Line)" in tab "Clip Object"
2. Place Landmarks at the model border, last point close to the first one
3. Press button "Trim" and confirm the correct trimming
4. Select and add socket in tab "Model Base"
5. Save as new finding

Module user interface



Links

	Handout 3D Scan Import
	https://youtu.be/TAAGkmFIKD4
	https://youtu.be/Q9R9V5QgEoE
	Import Adjust Scan

From:

<http://www.onyxwiki.net/> - [OnyxCeph³™ Wiki]

Permanent link:

http://www.onyxwiki.net/doku.php?id=en:add_image_3d&rev=1586879142

Last update: **2020/04/14 17:45**

