

Module Bonding Jigs 3D




Program version	BASIC	2D PRO	3D PRO	LAB	OMS
Availability	✗	✗			

In module Bonding Jigs 3D module, occlusal transfer jigs can be constructed to transference the virtually planned bracket positions using metallic standard bracket adapters.

First Steps

1. Open Module with bracket-attached Model
2. Select Crowns
3. Assign adapter and Jig parameters to crowns
4. Construct jigs individually or connected
5. Arrange jigs for manufacturing / export
6. Save planning as project and/or as finding

Module user interface

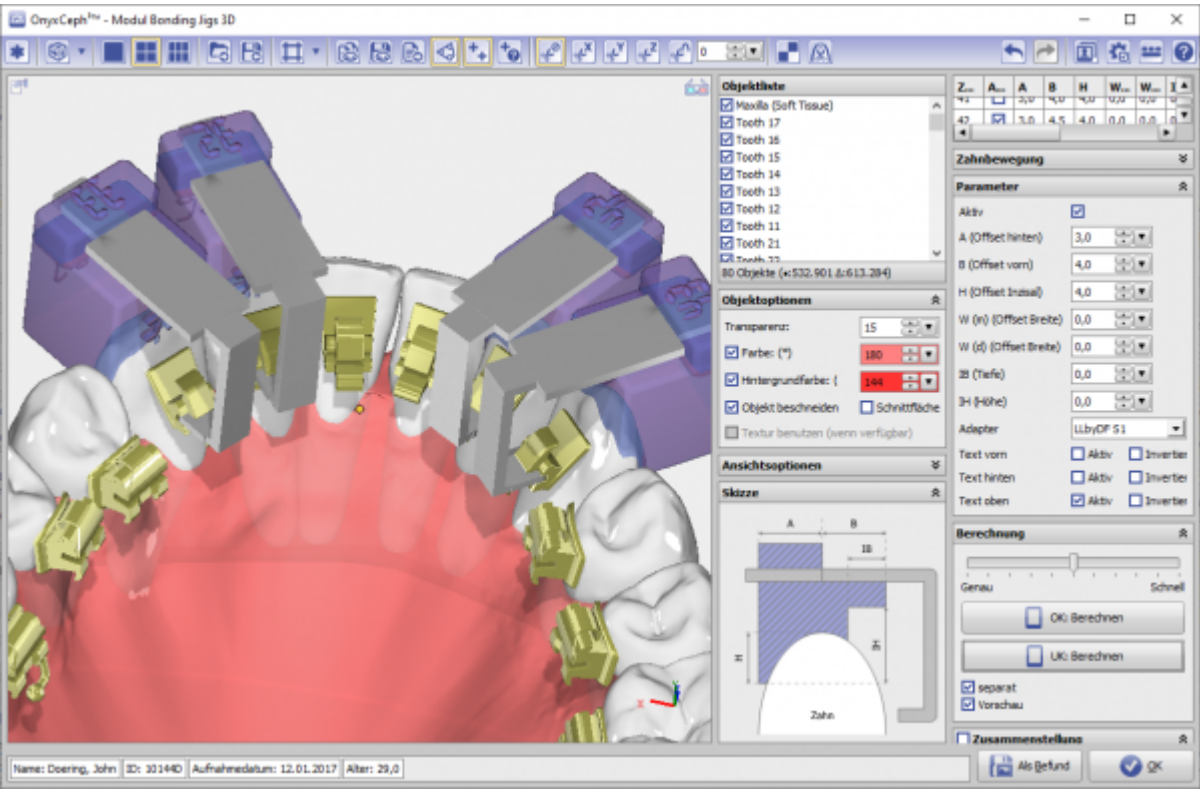

Links

[pdf.png_32](#)



From:

<http://www.onyxwiki.net/> - [OnyxCeph³™ Wiki]

Permanent link:

http://www.onyxwiki.net/doku.php?id=en:bonding_jigs_3d&rev=1613249945

Last update: **2021/02/13 21:59**

