


System Options - Tab |Third Party Tools|

On system options tab |Third Party Tools|, communication options can be configured.

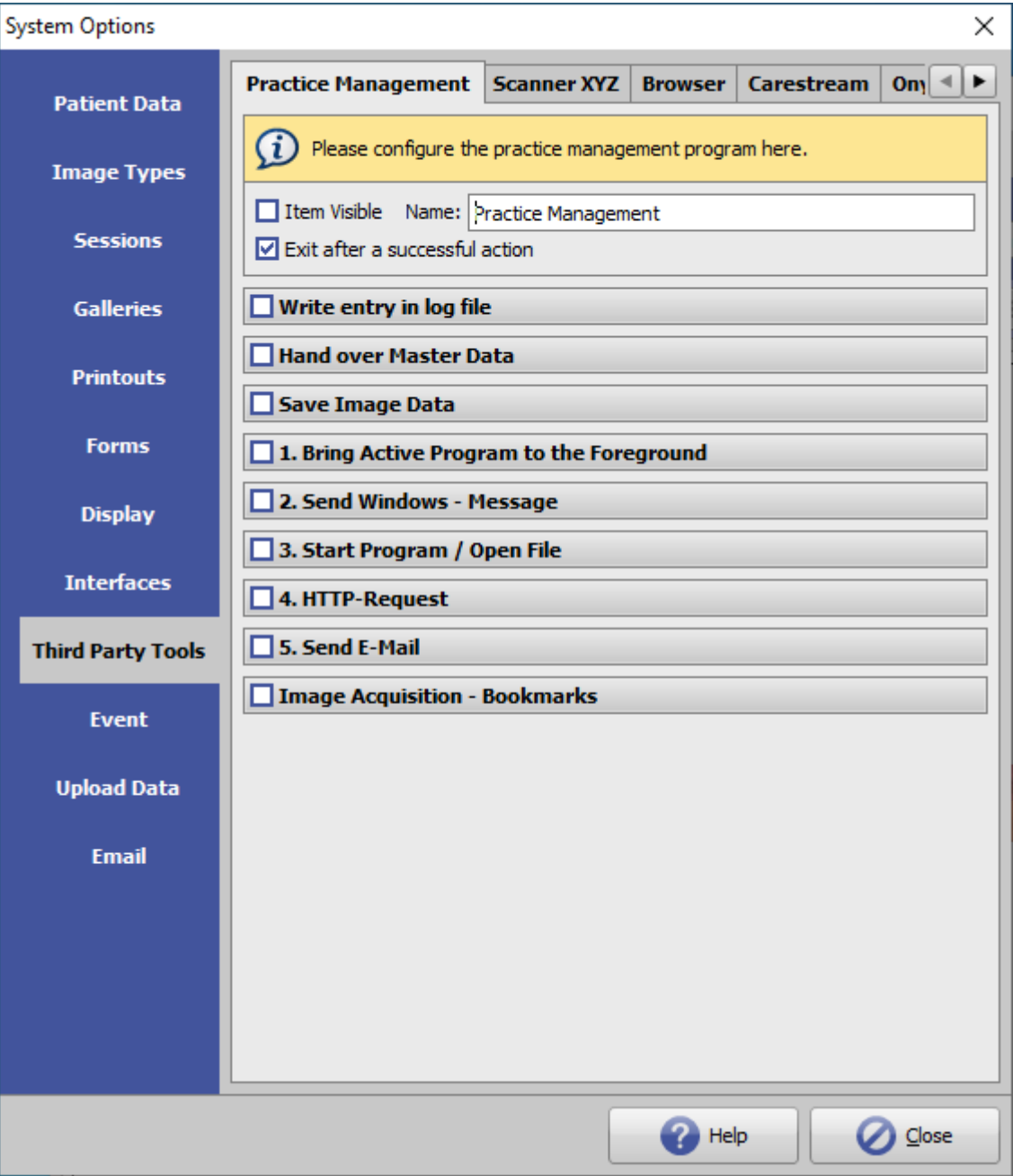

[Click on tab to switch to related topic]

Tools

A total of 5 third-party programs can be configured, which can be called via main menu [Go to]. Each call causes the optional reactions defined below.
If option [Image Acquisition - Bookmarks] is configured, the call is made in the bookmark.menu in

module Add image.

Options

Write entry in log file
Hand over Master Data
Save Image Data
1. Bing Active Program to the Foreground
2. Send Windows-Message
3. Start Program / Open File
4. HTTP Request
5. Send E-Mail
Image Acquisition - Bookmarks

If the option [] **Exit after a successful action** is not active, all active actions are executed.

Tab 1 is reserved for running the practice management software. Visibility and name of the related menu item can be set individually.

Example

 [Configuring the Windows™ webcam as an image source](#)

From:

<https://onyxwiki.net/> - [OnyxCeph³™ Wiki]

Permanent link:

https://onyxwiki.net/doku.php?id=en:tab_thirdpartytools

Last update:

2023/08/25 12:46

