


Interface Configuration Digital Xray: Cliniview

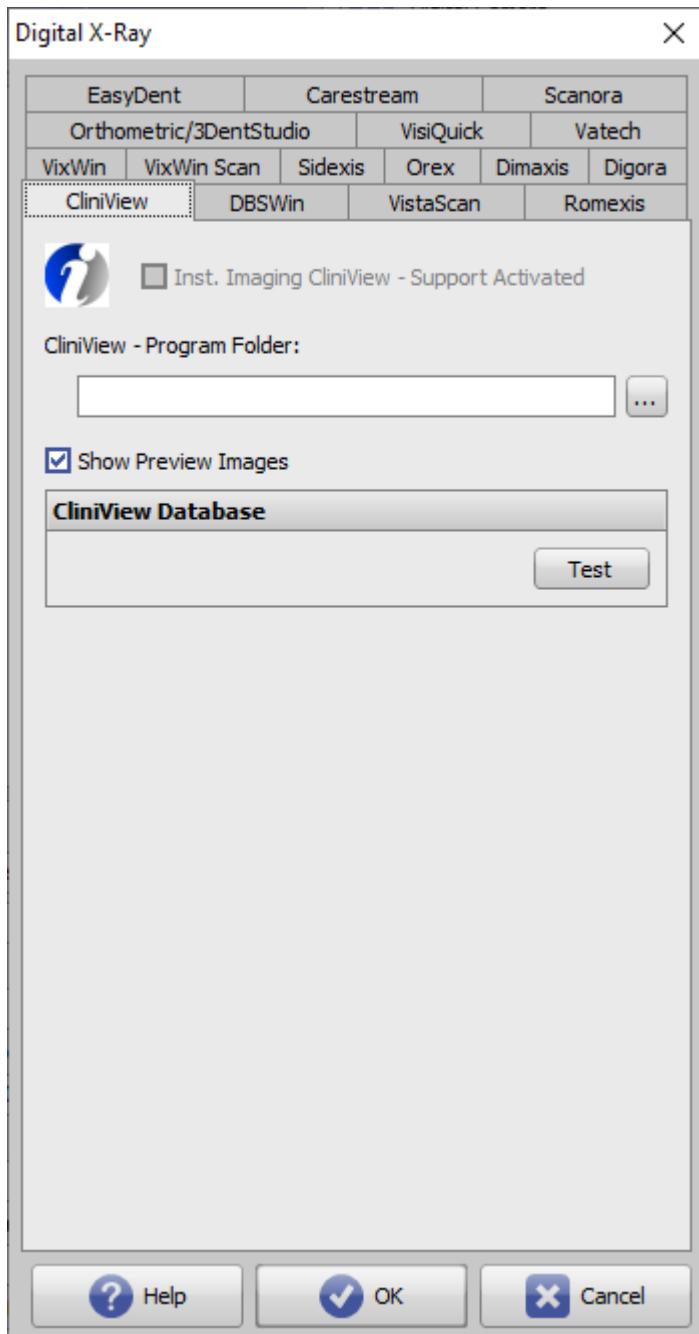

| Parameter | Setting |
|-----------------------|---|
| [Support Activated] | Activates import icon button in module Add Image 2D |
| Program Folder [...] | Browse for local application location ..\Cliniview.exe |
| [Show Preview Images] | Activate thumbnail preview in source preview panel in module Add Image 2D |
| [Test] | Testing database access |

From:

<http://onyxwiki.net/> - [OnyxCeph^{3™} Wiki]

Permanent link:

<http://onyxwiki.net/doku.php?id=en:cliniview>

Last update: **2022/06/27 07:32**

