


# Module FA\_Bonding

## First Steps

1. Place bracket point from list onto the correct tooth, model rotation determines up vector
2. Select bracket for point or load bracket set
3. Adjust bracket using the 3D navigators or the line-up tool
4. Harmonize crown distance
5. Save as project and new finding

## Module user interface

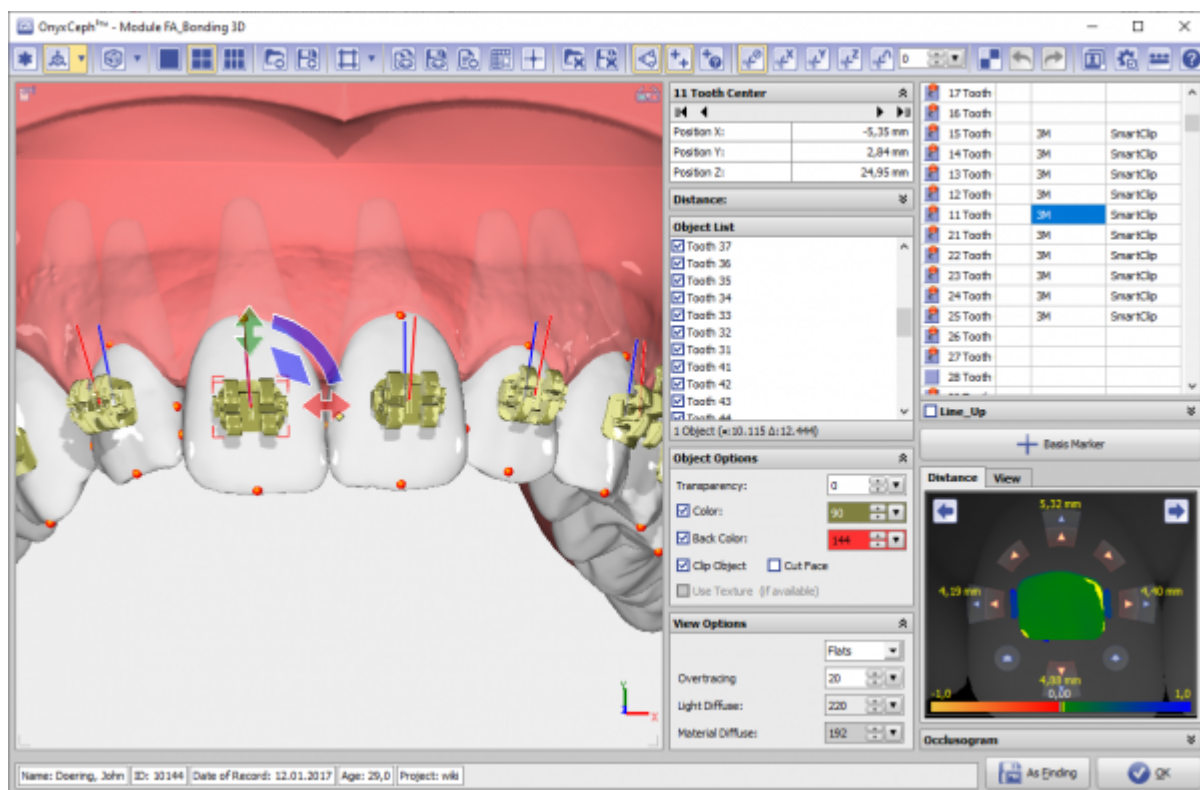

## Links



[https://image-instruments.de/Onyx3/Handout\\_IDB\\_FA\\_BT\\_DE.pdf](https://image-instruments.de/Onyx3/Handout_IDB_FA_BT_DE.pdf)



<https://youtu.be/bgCUllm2YKw>



[https://youtu.be/uulEFoQ\\_IV0](https://youtu.be/uulEFoQ_IV0)

From:

<http://www.onyxwiki.net/> - [OnyxCeph<sup>3</sup>™ Wiki]

Permanent link:

[http://www.onyxwiki.net/doku.php?id=en:fa\\_bonding&rev=1584709251](http://www.onyxwiki.net/doku.php?id=en:fa_bonding&rev=1584709251)

Last update: **2020/03/20 14:00**

