


# Interface Configuration Digital Xray: Romexis

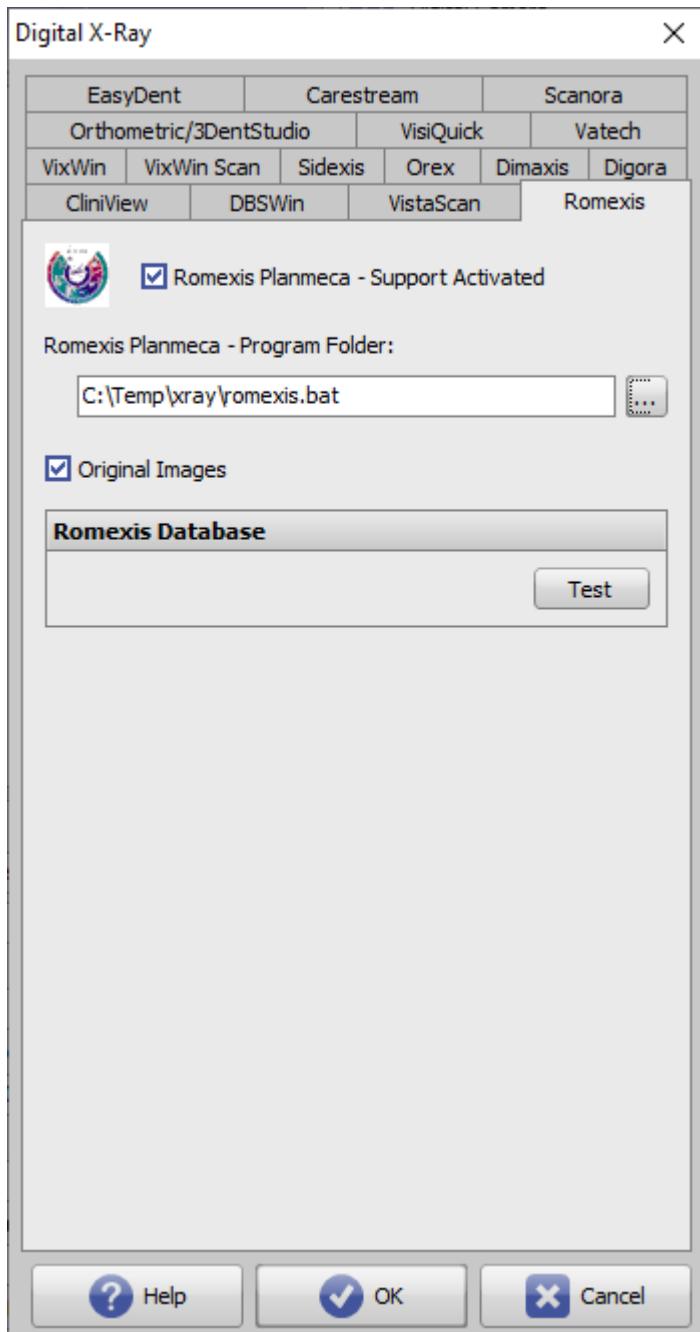

Parameter	Setting
[Support Activated]	Activates import icon button in module <a href="#">Add Image 2D</a>
Programmpfad [...]	Lokaler Pfad zur ausführbaren Datei ..\Romexis.bat
[Originalbilder]	Original- bzw. bearbeitete Bilder laden
[Test]	Test Datenbankzugriff

**Note:**

On the Romexis side, VDDS must be activated for the functionality of the interface and, if necessary, a special communication port (please take from the Romexis interface documentation) must be released.

## Known Problems

[DIDAPL.dll error](#)

From:  
<http://onyxwiki.net/> - [OnyxCeph<sup>3™</sup> Wiki]



Permanent link:  
<http://onyxwiki.net/doku.php?id=en:romexis>

Last update: **2023/10/09 11:43**