

Module Bonding Jigs 3D



<WRAP group> <WRAP threequarters column>


Program version	BASIC	2D PRO	3D PRO	LAB	OMS
Availability	×	×			

In module Bonding Jigs 3D module, occlusal transfer jigs can be constructed to transference the virtually planned bracket positions using metallic standard bracket adapters.

First Steps

1. Open Module with bracket-attached Model
2. Select Crowns
3. Assign adapter and Jig parameters to crowns
4. Construct Jigs individually or connected
5. Arrange Jigs for manufacturing / export
6. Save planning as project and/or as finding

Module user interface

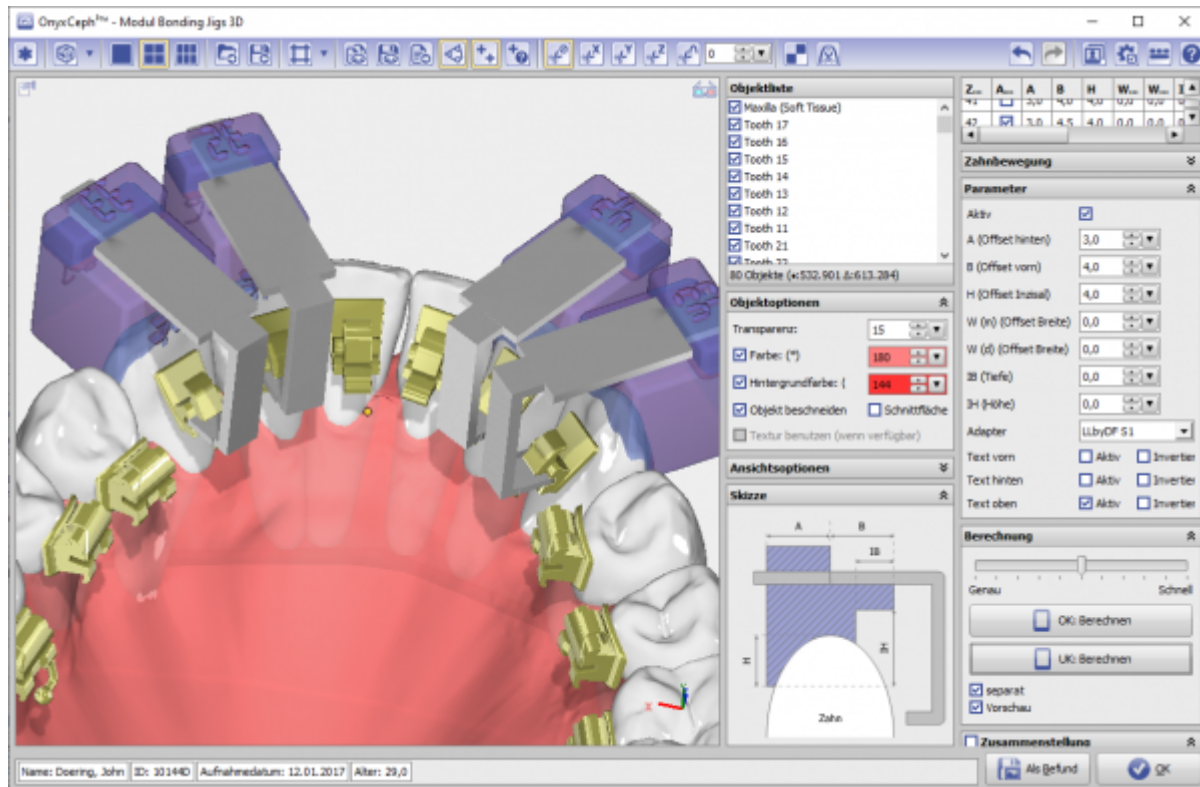

Links



From:
<https://www.onyxwiki.net/> - [OnyxCeph³™ Wiki]

Permanent link:
https://www.onyxwiki.net/doku.php?id=en:bonding_jigs_3d&rev=1644236344

Last update: **2022/02/07 13:19**

